

VFC400 Download of Data

Downloading data from your the VFC400

Control Solutions, Inc.
35851 Industrial Way, Suite D
St. Helens, OR 97051
503-410-5996
www.vfcdataloggers.com



Control Solutions, Inc.
888-311-0636

Downloading data

Buttons Actions Overview



Review/Mark button

Press to enter REVIEW mode and then press again to advance through the min/max temperatures. **This action also places an inspection mark in the log if RECORDING is active. This validates your twice daily temperature logs.**

Start/Clear/Stop button

Press to start recording
Press to Stop recording
Press to exit day summary



Downloading data

- The first thing we want to do before we download data from the VFC400 is stop the data logger from logging temperatures.



Downloading data

Stopping the data logger




When it is time to download the data from the data logger we suggest stopping the logger. To accomplish this you hold down the Start/Clear/Stop button for 1-2 seconds. The word STOPPED will appear on the screen as shown on the top left picture. The REC symbol will disappear after 1-2 seconds and you take your finger off the button quickly.

The logger will stop and the screen will only say STOPPED and show the number of days collected.



Downlaoding data

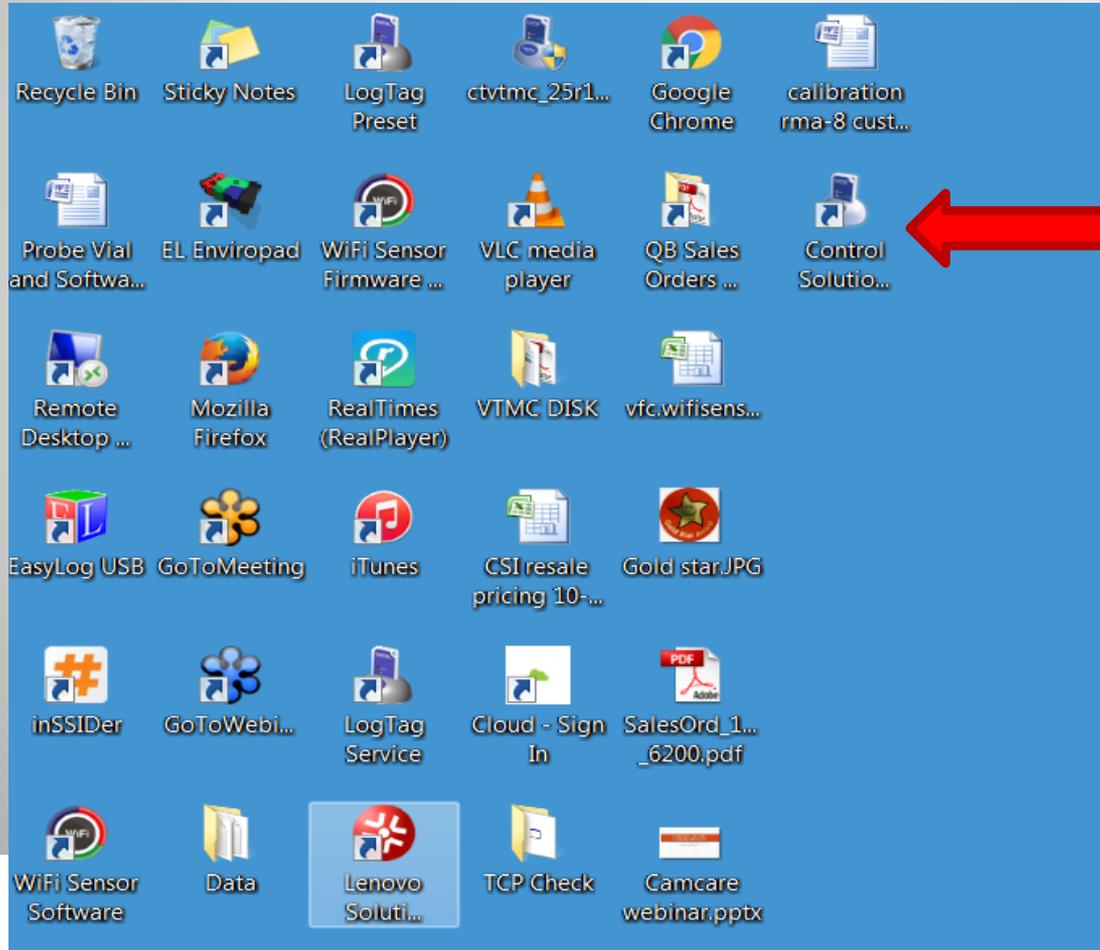

At your PC
open up the
Control
Solutions
VTMC
software by
double clicking
on the Icon.



VFC400 Data Logger Configuration

- Plug your Docking Station into the USB port of your PC
- Insert data logger into the Docking Station firmly




← Docking Station




Downloading data

- After you insert the logger into the docking station you will see the following message:

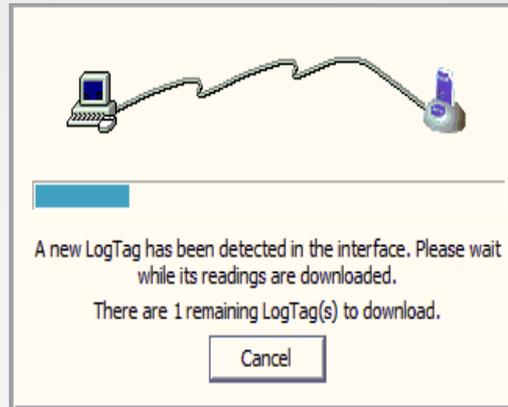

A few seconds later the chart/data will pop-up on the screen and then the following message:


Click OK in the box and remove the data logger from the docking station.

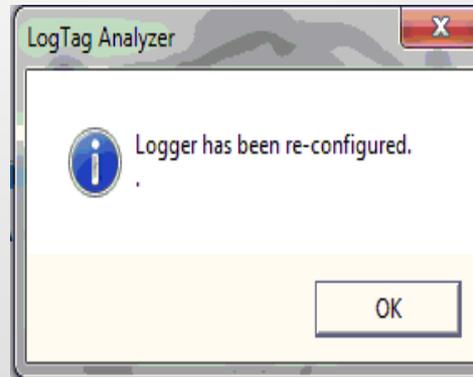


Downloading data

Starting The Logger

Recorder in *READY* mode

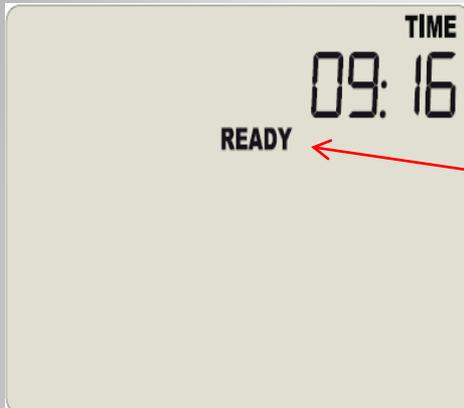

The logger is now ready to start.
In Ready mode the recorder displays the time and the READY icon.

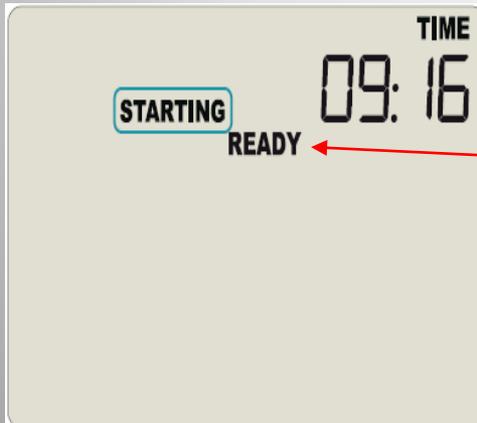


To Start: Press and hold the Start/Clear/Stop button. First, STARTING is displayed in addition to READY. Then READY disappears.



Downlaoding data

Starting The Logger

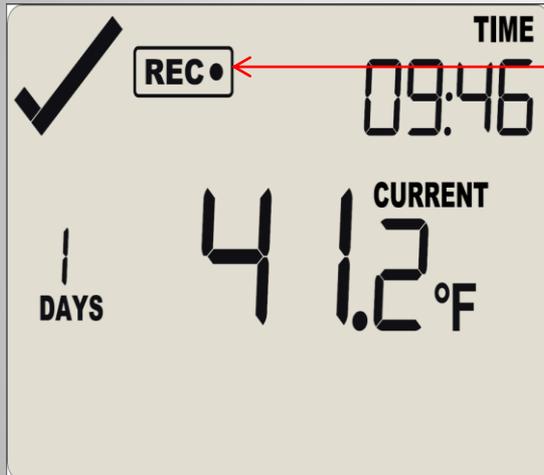

Once READY disappears, immediately release the button. STARTING will also disappear, and the symbol REC• will be shown. The VFC 400 is now recording temperature data.

You can see that there are no alarms because there is a ✓ mark in the upper left corner and the audible alarm is not beeping.





Control Solutions, Inc.
www.vfcdataloggers.com
888-311-0636

Thank you for your business